

BUTTERFLY INSTRUCTIONS



01 OPEN

Remove your butterfly wrap and accessory straps from their original packaging.



02 PLACE

Place the butterfly wrap over the therapy site according to your doctor's instructions.



03 SECURE

Secure the butterfly wrap as needed, using the accessory straps.



04 COMFORT

Once fitted, the wrap should be held snugly in place with no creases or folds.

PRODUCT AND WARRANTY INFORMATION

<https://www.nicerecovery.com/instructions-for-use>

**THIS
IS NICE**

For more information:
Call a customer representative
at (888) 815-9907 or email
info@nicerecovery.com

Location:
Boulder, CO 80301



Visit:
nicerecovery.com



Outside the United States: Please contact your local distributor.



6-1016-01 REV. A

Consult a healthcare professional before using this product




- Connect the power cord to the external power supply.
- Plug the circular connector into the back of the Cooling Unit. You will hear a click when properly connected.
- Plug the power cord into the wall outlet.



- Connect the hose to the control unit. You will hear a click when properly connected.



- Place the wrap on the appropriate body part. Secure with hook fastener strips.
- Connect the hose to the wrap. You will hear a click when properly connected.

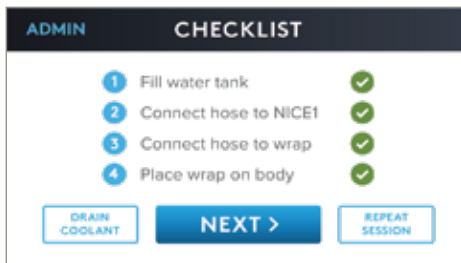

STEP 1

Turn on the power switch on the back of the device. The touch screen will illuminate and display the "Nice" logo. Press Start. Review the checklist on the touch screen. Confirm and press NEXT.




STEP 2

Select desired COLD LEVEL and press NEXT. (level 5 is the coldest setting).


Cold settings are as follows:

Level 1 = 58F (14.4C) Level 4 = 46F (7.7C)
Level 2 = 54F (12.4C) Level 5 = 42F (5.5C)
Level 3 = 50F (10C)

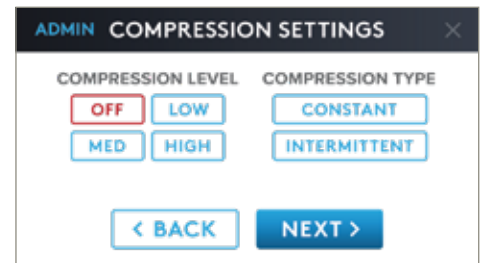

STEP 3

Select desired COMPRESSION LEVEL and COMPRESSION TYPE and press NEXT.

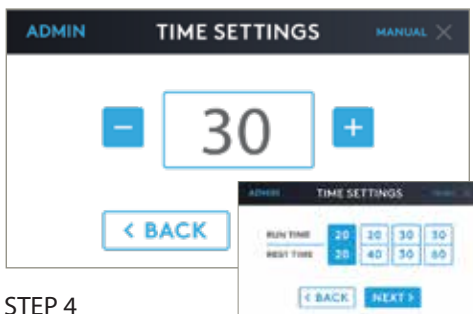


STEP 4

Choose a MANUAL time setting of 5 – 40 minutes for a single therapy session OR choose a PRESET PROGRAM for multiple sessions and press NEXT. Note: Preset timing programs are on a 10-cycle memory.



STEP 5

Review and confirm your settings by pressing START or press BACK to change the settings.

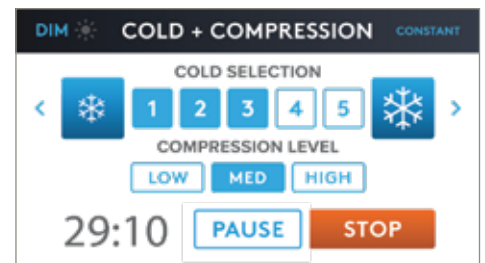


STEP 6

During your therapy session you can change the cold or compression settings simply by pressing 1 – 5 or HIGH, MED or LOW. To Dim, Press Dim switch, to brighten — touch screen.



Caution: Consult Accompanying Documents.

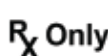

Caution: Federal law restricts this device to sale by or on the order of a physician.



Medical Device



Niche Recovery Systems LLC
2205 Central Ave, Suite A
Boulder, CO 80301



Date of manufacture



For more information call a customer service representative at (888)815-9907 or email info@nicerecovery.com



Non-sterile



Do not use if package is damaged and consult instructions for use



Does not contain natural rubber latex



Single patient multiple use



MDSS GmbH
Schiffgraben 41
30175 Hannover, Germany